

The New Movement In Orthodontics™



NeoMaestro™ Premium Bracket System

Bracket System
A Family of Brackets that
All Work in Unison



# Strategy. Research. Evolution. NeoMaestro Premium Bracket Systems

In 2009 our vision began to develop the ultimate twin-bracket system. One where the brackets were interchangeable between aesthetic and metal for convenient and cost-efficient bracketing. This unique system would provide clinicians optimal benefits and conveniences, while delivering patient comfort and successful treatment outcomes. The Company's Research and Development team spent laborious time analyzing twin brackets currently available in the market, as well as soliciting feedback from clinicians. The result is a modern day success: *NeoMaestro* Family of Brackets all designed for interchangeability.

The three twin-bracket lines offer uniform design features that are unique to only these systems. They allow a doctor to utilize aesthetic brackets in the anterior and metal brackets in the posterior without comprising patient treatment.

- **Comparable bracket sizes** Optimal bracket interchangeability
- Compatible In/Out No wire compensating bends required between bracket transitions
- Same bracket base compound contours Consistent adaptation to the teeth
- Equivalent bond strength Consistent and reliable debonding between bracket systems
- Three twin-bracket material options Ability to cater to patient/doctor preferences and save money
- Harmonious Maestro Buccal Tubes Complete system approach
- **Consistent prescriptions** Ability to meet clinician preferences on torque & tip, and can mix and match to meet clinician preferences on Rx values
- **Uniform color-coding –** Easy bracket identification

#### COMPARABLE BRACKET SIZES

The family of brackets were designed and developed with a common size in mind to give you the same feel and handling as you work with each bracket.



Maestro® Bracket
Premium Stainless Steel



**NeoLucent** <sup>®</sup> **Plus**<sup>™</sup> **Bracket** Polycrystalline



NeoCrystal™ Bracket
Monocrystalline

## **COMPATIBLE IN/OUT**

The complete *NeoMaestro* Family of Brackets has compatible in/out's that allow you to easily mix our aesthetic *NeoLucent Plus* Ceramic Brackets or *NeoCrystal Sapphire* Ceramic Brackets with our *Maestro* Metal-Twin Brackets. Due to the exact tolerances between the brackets, there is no need for wire bends as you transition from the anterior aesthetic brackets to the metal brackets in the posterior.



## SAME BRACKET BASE COMPOUND/CONTOURS

This ensures that you experience the same optimized fit to tooth morphology no matter if you are working with a *Maestro, NeoLucent Plus,* or *NeoCrystal* Bracket.








Maestro Bracket


NeoLucent Plus Bracket

NeoCrystal Bracket

## **EQUIVALENT BOND STRENGTH**

Whether you are working with the *Maestro* Bracket's 80-gauge mesh pad, *NeoLucent Plus* or *NeoCrystal* Bracket's propriety *Crunch COAT™* Base, you will have the same bonding experience. Our propriety *Crunch COAT* Base behaves similar to our 80-gauge mesh pad. In addition, both bases provide a easy surface area for bond strength.

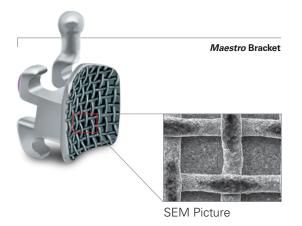

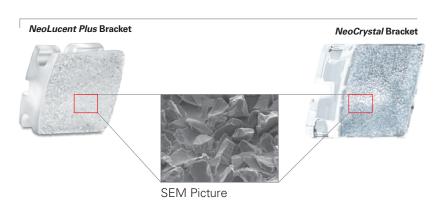

## THREE TWIN-BRACKET MATERIAL OPTIONS

The three style of brackets that compose the *NeoMaestro* Bracket System are built upon a true twin-bracket design with torque in base.



**Maestro Bracket** Premium Stainless Steel



**NeoLucent Plus Bracket** Polycrystalline



**NeoCrystal** Bracket Monocrystalline

# HARMONIOUS MAESTRO MOLAR BUCCAL TUBES

*Maestro* Buccal Tubes complete the *NeoMaestro* Family of Brackets. They have a large funnel entrance, side-grip areas, and color-coding to make for easier placement, position, and identification. All three bracket systems can be utilized with our bondable or weldable 1st and 2nd molars for optimized inventory.



Maestro 1st Molar Buccal Tube



Maestro 2<sup>nd</sup> Molar Buccal Tube


## **CONSISTENT PRESCRIPTIONS**

The three bracket systems are offered in the same prescriptions so that brackets can easily be mixed and matched, while still resulting in the same optimal treatment outcome that you are accustomed to receiving.

	CENTRALS		LATERALS			CUSPIDS			1ST BICUSPIDS			2ND BICUSPIDS			
	TRQ.	ANG.	ANTI- ROT.	TRQ.	ANG.	ANTI- ROT.	TRQ.	ANG.	ANTI- ROT.	TRQ.	ANG.	ANTI- ROT.	TRQ.	ANG.	ANTI- ROT.
MAXILLARY															
McLaughlin, Bennett, and Trevisi	+17°	+4°	0°	+10°	+8°	0°	-7°	+8°	0°	-7°	0°	0°	-7°	0°	0°
McLaughlin, Bennett, and Trevisi (0 torque on 3's)	+17°	+4°	0°	+10°	+8°	0°	0°	+8°	0°	-7°	0°	0°	-7°	0°	0°
Roth	+12°	+5°	0°	+8°	+9°	0°	-2°	+13°	4°M	-7°	0°	2°D	-7°	0°	2°D
MANDIBULAR															
McLaughlin, Bennett, and Trevisi	-6°	0°	0°	-6°	0°	0°	-6°	+3°	0°	-12°	2°	0°	-17°	2°	0°
McLaughlin, Bennett, and Trevisi (0 torque on 3's)	-6°	0°	0°	-6°	0°	0°	0°	+3°	0°	-12°	2°	0°	-17°	2°	0°
Roth	-1°	0°	0°	-1°	0°	0°	-11°	+7°	2°M	-17°	0°	4°D	-22°	0°	4°D

### UNIFORM COLOR-CODE IDENTIFICATION

Easily identify brackets between the *NeoMaestro* Bracket Systems, as similar color-code ID markings are used. Looking at a black ID dot on the distal gingival tie-wing of either a *Maestro, NeoLucent Plus*, or *NeoCrystal* Bracket, you will confidently know it is an upper right central.

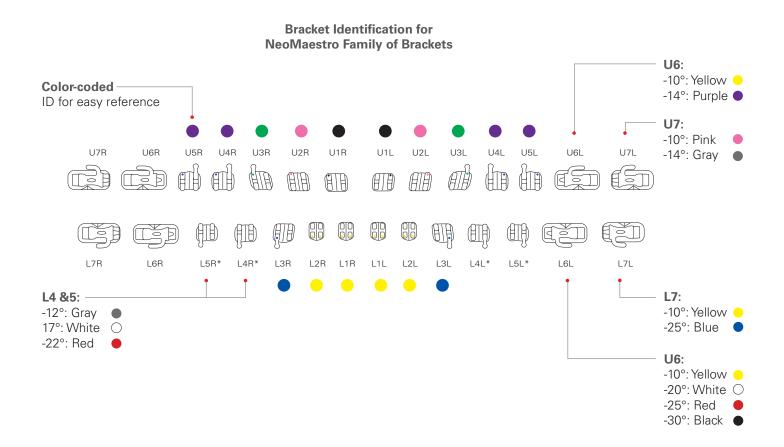

# OPTIONS FOR METAL/AESTHETIC BRACKET KITS

To help simplify your ordering process we have developed 5x5 kits that are composed of both aesthetic and metal brackets. You will receive aesthetic brackets, based on your choice of either *NeoCrystal* or *NeoLucent Plus* Ceramic Brackets, for the 3x3 section, and *Maestro* Brackets for the bicuspids.

	NEOCRYSTAL/MAI	ESTRO KITS	NEOLUCENT PLUS/MAESTRO KITS			
Rx	.018	.022	.018	.022		
McLaughlin, Bennett, and Trevisi	714-296MC-NM	714-396MC-NM	704-296MC-NM	704-396MC-NM		
McLaughlin, Bennett, and Trevisi (0 torque on 3's)	714-296MC-NM	714-196MC-NM	-	-		
Roth	714-096-NM	714-396-NM	704-296-NM	704-396-NM		



All kits U/L 5x5 (3x3 aesthetic brackets with metal brackets on 4's and 5's. Hooks on 3's, 4's, and 5's)\*

The mission of *Henry Schein Orthodontics* is to provide state-of-the-art orthodontic products and innovative clinical solutions that enable our customers to offer exceptional patient care, while expanding the scope and profitability of their practices.



For more information on our products and educational offerings, please contact us: HenryScheinOrtho.com
Customer Service: 888.851.0533
In the U.S.: 800.547.2000

Outside the U.S.: +(1) 760 448 8600 Canada: CERUM 800.661.9567 To email an order: USASales@HenryScheinOrtho.com To email an order outside the U.S.: IntlOrders@HenryScheinOrtho.com







The New Movement In Orthodontics The

Ortho Organizers, Inc. 1822 Aston Ave. Carlsbad, CA 92008-7306 USA