





Racegel

Full power for gingival management



RACEGEL IS THE ONLY THERMO-GELIFIABLE GEL CREATING A CLEAN AND DRY ENVIRONMENT FOR ANY PROCEDURE WHERE GINGIVAL TISSUE MANAGEMENT IS NECESSARY

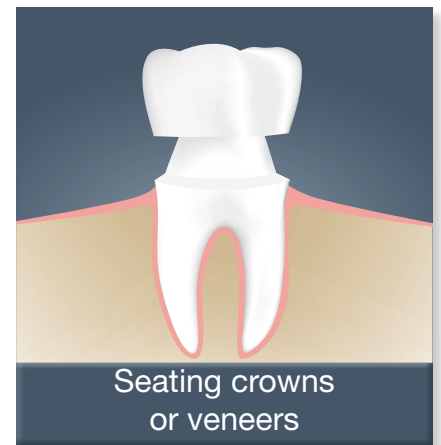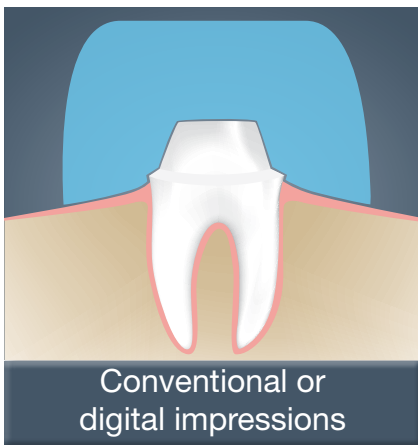


Integrated Astringent Effect

- **Soft tissue contraction** thanks to 25% aluminum chloride
- **Optimal control** of bleeding without staining of the gingiva
- **Easy access** to all subgingival margins
- Specific hemostatic solutions **not required**
- Creation of a **clean and dry** environment



Atraumatic Gingival Preparation

- **From liquid to gel** : Unique thermo-gelifiable viscosity:
 - At room temperature (20°C), Racegel is liquid in the syringe
 - Racegel viscosity increases with temperature
 - In contact with oral tissues (35°C) Racegel viscosity increases and immediately becomes a gel

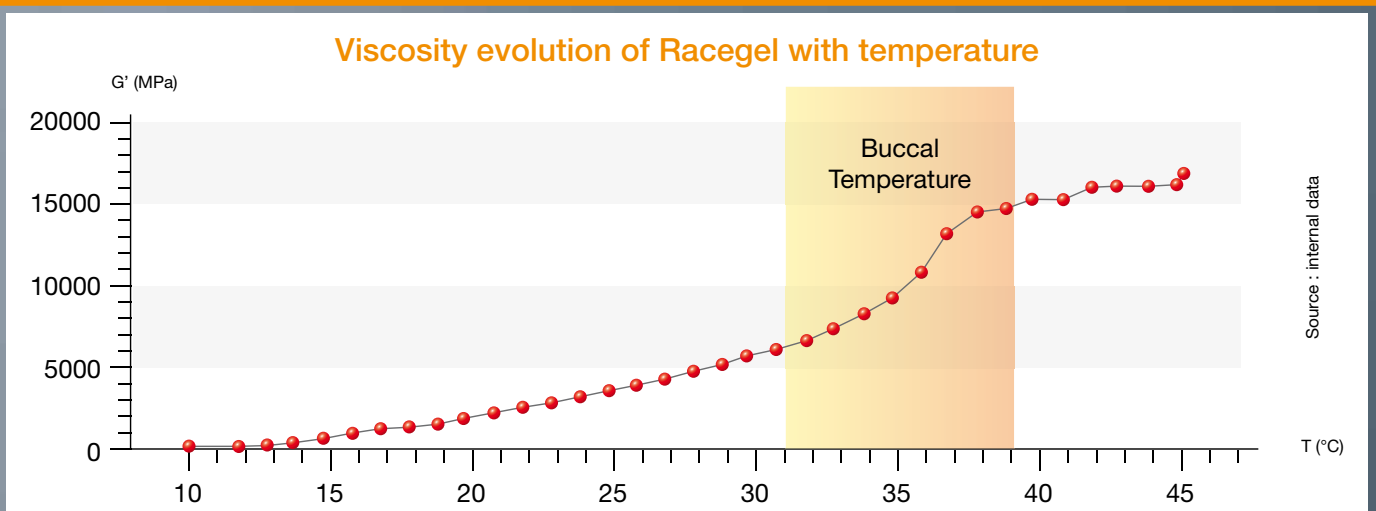

Easy, Fast and Precise Placement

- **Ready-to-use syringe** with pre-bent needle tips
- **Easy insertion** of the gel into the sulcus
- The gel does **not run or drip** during placement
- **Easy visualization of the gel** thanks to its orange color



Optional use with retraction cords in case of a deep subgingival preparation

Fast and Secure Rinsing

- **Reversible** thermo-gelifiable viscosity:
 - Racegel viscosity decreases when rinsed with cold water
 - Absence of rebleeding when Racegel is rinsed away



Minimal displacement yet high quality impressions




Minimal yet sufficient gingival displacement for accurate record of the cervical margins for high quality impressions every time

Racegel

the efficient way to manage
all your gingival preparations

Racegel vs. Retraction Pastes	Racegel vs. Hemostatic Solutions	Racegel vs. Retraction Cords
Does not require additional hemostatic solution	Stays in place, does not run or drip	Is easier and faster to put in place
Ensures high quality impressions with minimal gingival displacement	Does not cause any gingival staining	Does not cause rebleeding during removal
Spars the need for dedicated equipment	Ensures high quality impressions with a minimal gingival displacement	Ensures high quality impressions with a minimal gingival displacement
Has a raspberry flavor for patient comfort	Offers a convenient ready-to-use delivery system	Maintains the integrity of the junctional epithelium

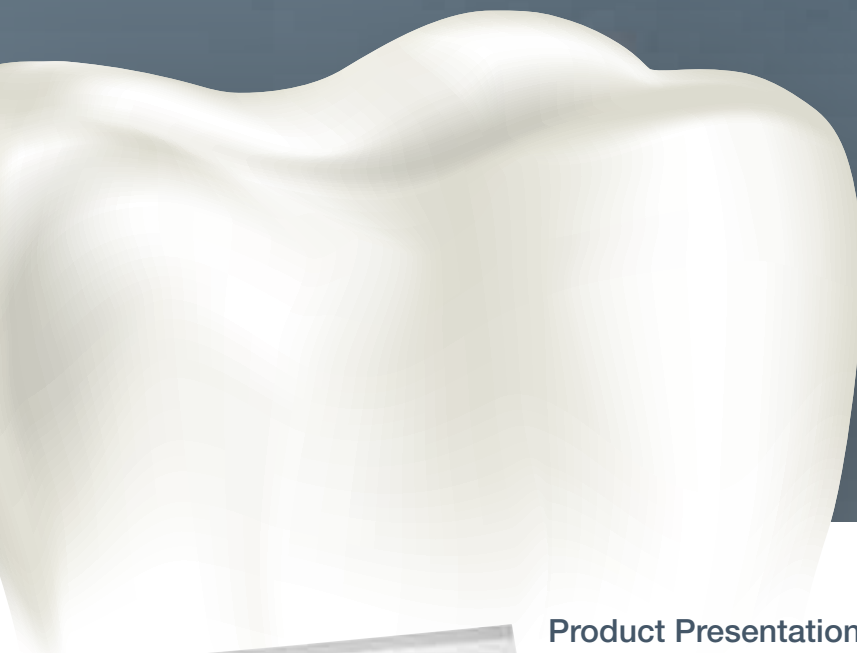

Product Presentation



Available in:

- Box of three syringes and 30 pre-bent tips



The Dental Pharmaceutical Company

INNOVATIVE, SAFE AND EFFECTIVE SOLUTIONS FOR DENTISTRY WORLDWIDE